# POLARISOL Pro

# POLARISOL Mobilo

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Welcome to the world of POLARISOL, congratulations on your decision. We are confident that POLARISOL products accompany you in your everyday life, while preserving the health and beauty of you and your loved ones.

Please read carefully the operating instructions, handling and security specifications before use! Keep the operating instructions in a safe place for future reference! Do not give the device to anyone without the operating instructions!

## POLARISOL – polarised light treatment every day

POLARISOL is a unique regenerative appliance with a special and innovative optical unit and a large light -emitting surface, producing higher quantities of concentrated and invigorating polarised light, that is difficult to access in nature.

POLARISOL is a modern and improved version of a 30-year-old technological solution and has been developed for general regeneration and feel-good purposes.

POLARISOL is also suitable for home and professional use and has been created by the rapid development of modern technology.

It does not contain ultraviolet (UV) rays and the degree of polarisation of its light is more than 99%.

POLARISOL is not a medical device and it is not intended for preventing, diagnosing, monitoring, treating a disease or alleviating the symptoms of a disease. It is not suitable for the diagnosis, monitoring, treatment or alleviation of injury or handicap. It does not investigate, replace, compensate or modify the anatomy or a physiological process. POLARISOL has been exclusively developed for regeneration, well-being and cosmetic purposes.

The quality technology and innovative design combined with the ease of use brings the result that POLARISOL, thanks to polarised light, is able to enhance the body's energy supply, optimizes its operation, protects against harmful external influences, balances processes back to steady state and increases cell vitality.

In addition it is ideally used for cosmetics treatments, as it is effective in the rejuvenation of the skin and in the fight against wrinkles and cellulite.

The POLARISOL stand applicable to the appliance allows an easy setting in various positions, and ensures a convenient light treatment of any part of the body.

When using POLARISOL, keep a distance of 10 cm from the surface treated! The recommended duration of a treatment is 15 to 20 minutes, up to several times a day (2-3 occasions).

#### **It is recommended for:**

beauty, wellness, feel-good, energizing, regeneration, revitalizing skin and mood-enhancing purposes.

#### **The characteristics of POLARISOL polarised light:**

- the light emitted by POLARISOL devices is  $\geq 99\%$  polarised
- polarised light waves oscillate on parallel planes
- POLARISOL light is polychromatic light, containing several range of wavelengths, including visible light and a part of the lower infrared spectrum
- at low temperature it penetrates into the skin without causing any damage in skin tissue

POLARISOL is not a medical device and it is not intended for preventing, diagnosing, monitoring, treating a disease or alleviating the symptoms of a disease. It is not suitable for the diagnosis, monitoring, treatment or alleviation of injury or handicap. It does not investigate, replace, compensate or modify the anatomy or a physiological process. POLARISOL has been exclusively developed for regeneration, well-being and cosmetic purposes.

- it is similar to natural sunlight, but does not contain harmful ultra violet (UV) radiation
- unlike the laser light POLARISOL is an incoherent light (out of phase) and as in the case of natural sunlight, the light waves are not synchronized
- a combination of innovative Hungarian patent and Japanese technology
- versatile
- characterised by a short treatment time
- easy to use
- light penetrates optimally the tissues
- it's minimal thermal effect does not damage the skin and the underlying tissues

#### **The beneficial effects of POLARISOL:**

- support of human defense system
- increase of general well-being
- relieve of stress and anxiety caused by lack of light
- improving the ability to concentrate
- increasing physical and mental performance
- regeneration
- support of muscle strenght and capacity
- cellulite treatment (stimulation of cell metabolism, blood and lymphatic circulation)
- skin and hair care (improvement of skin and scalp oxygenation)
- improvement of locomotor functions

#### **Important! Do not use POLARISOL in case of any of the following symptoms:**

- photoallergy
- porphyria (all types)
- inflammation of the eye
- retinal diseases
- photosensitivity induced by drugs/herbs



POLARISOL is not a medical device and it is not intended for preventing, diagnosing, monitoring, treating a disease or alleviating the symptoms of a disease. It is not suitable for the diagnosis, monitoring, treatment or alleviation of injury or handicap. It does not investigate, replace, compensate or modify the anatomy or a physiological process. POLARISOL has been exclusively developed for regeneration, well-being and cosmetic purposes.

# Description of POLARISOL Pro

An ideal device both for family, home and institutional (eg.: therapeutic spa and wellness centers, thermal baths, beauty salons, sports clubs, hairdressing salons) use with a polarised light emitting optical surface of 402 cm<sup>2</sup>.

- ① Multi-layer film-based polarised light optical unit
- ② Display and control unit (Set time with +/- button to a maximum of 99 minutes)
- ③ Cable socket/connector
- ④ On/off switch
- ⑤ Backside ventilation slits

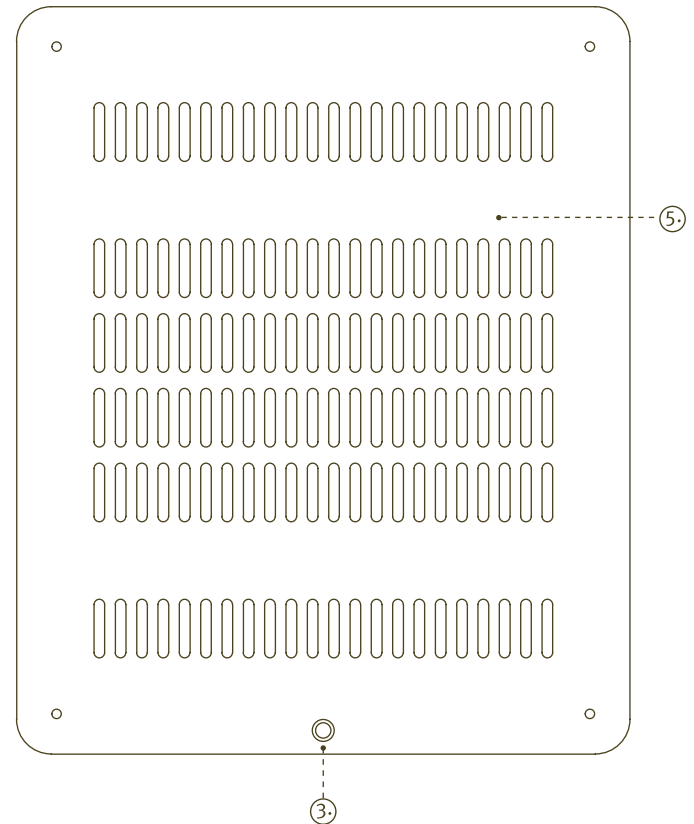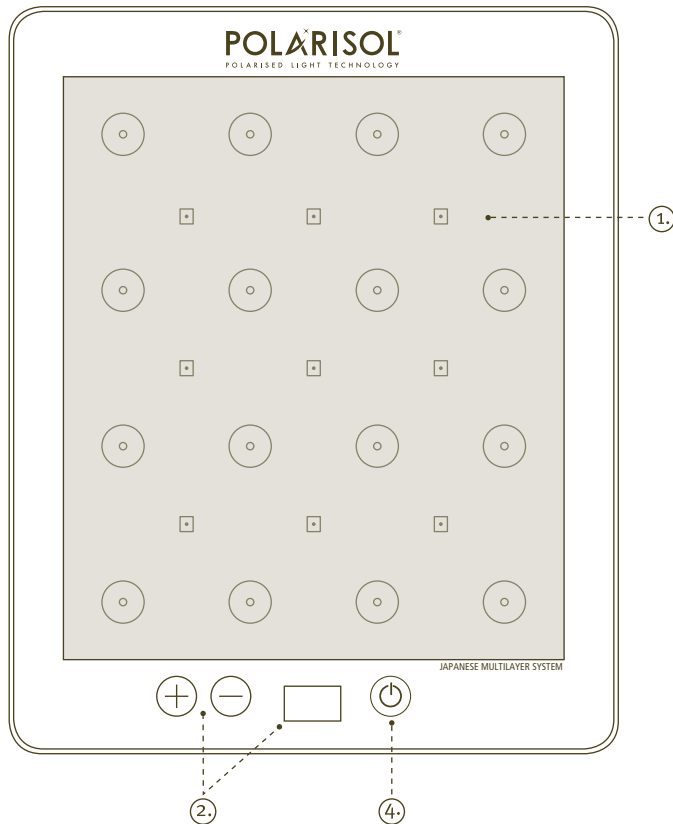POLARISOL is not a medical device and it is not intended for preventing, diagnosing, monitoring, treating a disease or alleviating the symptoms of a disease. It is not suitable for the diagnosis, monitoring, treatment or alleviation of injury or handicap. It does not investigate, replace, compensate or modify the anatomy or a physiological process. POLARISOL has been exclusively developed for regeneration, well-being and cosmetic purposes.

POLARISOL is not a medical device and it is not intended for preventing, diagnosing, monitoring, treating a disease or alleviating the symptoms of a disease. It is not suitable for the diagnosis, monitoring, treatment or alleviation of injury or handicap. It does not investigate, replace, compensate or modify the anatomy or a physiological process. POLARISOL has been exclusively developed for regeneration, well-being and cosmetic purposes.

## POLARISOL floor stand

In case you are using the appliance with a POLARISOL table/ floor stand, always place it on a dry, solid and horizontal surface.

When it is hold in hands, be sure to grasp it firmly with a dry (without water or other moisture) palm and fingers!  Protect the device from falling, scratching and hitting against to sharp or hard objects or surfaces.

The backside ventilation slints should be never blocked or covered by hands or by other materials because if the device ventilation is inadequate, then it shuts down automatically. In order to protect from dust and dirt, please store and transport the device and the cable cord in the travel bag designed especially for this purpose, that you will find in product packaging.

## Using POLARISOL Pro


①. Connect the device to the network. 

②. Clean and expose the body surface to be treated! Set the desired treatment time to the minute (the + button is to increase the – button is to decrease the treatment time), and then turn the device on by pressing the ⏻ button. The duration of the treatment can be stopped at any time prior to finishing the treatment by pressing the ⏸ button. 

POLARISOL is not a medical device and it is not intended for preventing, diagnosing, monitoring, treating a disease or alleviating the symptoms of a disease. It is not suitable for the diagnosis, monitoring, treatment or alleviation of injury or handicap. It does not investigate, replace, compensate or modify the anatomy or a physiological process. POLARISOL has been exclusively developed for regeneration, well-being and cosmetic purposes.

The unit will count backwards the treatment time, so you can check treatment.

③. Do not move during treatment! Polarized light should reach the treated body surface evenly from a distance of 10 cm at a right angle of 90°! In case of treating larger body areas treat them point-by-point! Treat the desired body surface for 15-20 minutes 2-3 times a day!

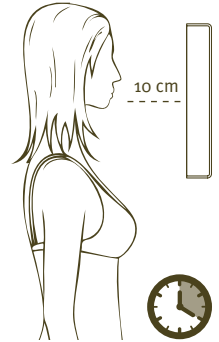The POLARISOL Pro stand allows POLARISOL to be used comfortably in any treatment positions.

④. At the end of the preset treatment time the device turns off automatically, or you can direct the light of the device to another body surface.


Always unplug the device after use! 

POLARISOL is not a medical device and it is not intended for preventing, diagnosing, monitoring, treating a disease or alleviating the symptoms of a disease. It is not suitable for the diagnosis, monitoring, treatment or alleviation of injury or handicap. It does not investigate, replace, compensate or modify the anatomy or a physiological process. POLARISOL has been exclusively developed for regeneration, well-being and cosmetic purposes.

# Description of POLARISOL Mobilo

The light emitting optical surface of POLARISOL Mobilo is 43 cm<sup>2</sup>. It is primarily recommended for use at home, but we also propose for beauty industry professionals as well.

Thanks to its ergonomic design and light weight, it is easy to use and carry for all ages.

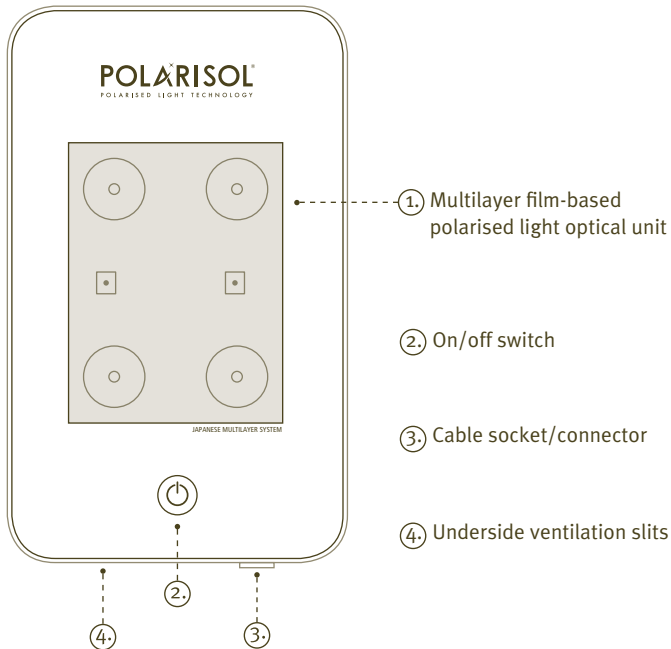POLARISOL is not a medical device and it is not intended for preventing, diagnosing, monitoring, treating a disease or alleviating the symptoms of a disease. It is not suitable for the diagnosis, monitoring, treatment or alleviation of injury or handicap. It does not investigate, replace, compensate or modify the anatomy or a physiological process. POLARISOL has been exclusively developed for regeneration, well-being and cosmetic purposes.

In case you are using the appliance with a table stand, always place it on a dry, solid and horizontal surface!

When it is hold in hands, be sure to grasp it firmly with a dry (without water or other moisture) palm and fingers!

Protect the device from falling, scratching and hitting against to sharp or hard objects or surfaces!

The underside ventilation slints should be never blocked or covered by hands or by other materials because if the device ventilation is inadequate, then it overheats and shuts down. In order to protect from scratches, dust and dirt, please store and transport the device and the cable cord in the travel bag designed especially for this purpose, that you will find in product packaging!

## Using POLARISOL Mobilo

①. Connect the device to the power outlet.

②. Clean and expose the body surface to be treated!  
The device will start and stop by pressing the  button.

For immediate re-start after switching off the device please wait a few seconds for switching on the LED bulbs!

③. Do not move during treatment!  
Polarized light should reach the treated body surface evenly from a distance of 10 cm at a right angle of 90°. In case of treating larger body areas move from point-by-point!

④. Treat the desired body surface for 15-20 minutes 2-3 times a day!

POLARISOL is not a medical device and it is not intended for preventing, diagnosing, monitoring, treating a disease or alleviating the symptoms of a disease. It is not suitable for the diagnosis, monitoring, treatment or alleviation of injury or handicap. It does not investigate, replace, compensate or modify the anatomy or a physiological process. POLARISOL has been exclusively developed for regeneration, well-being and cosmetic purposes.

At the end of the treatment you can finish it by by turning device off or direct the light of the device to another body surface.

- ⑤ Always unplug the device after use!

POLARISOL is not a medical device and it is not intended for preventing, diagnosing, monitoring, treating a disease or alleviating the symptoms of a disease. It is not suitable for the diagnosis, monitoring, treatment or alleviation of injury or handicap. It does not investigate, replace, compensate or modify the anatomy or a physiological process. POLARISOL has been exclusively developed for regeneration, well-being and cosmetic purposes.

## OPERATION/SETUP GUIDE

### Prior to operation

A 230V grounded socket/connector is required for the operation of the device

**THE DEVICE CAN BE USED ONLY WITH GROUNDED ELECTRIC NETWORK SYSTEM!**

The device in all cases should be placed in a way, that the electric connector (plug) is accessible at any time!

### Intended use

Each device is intended for indoor use only. Operate the unit only in non-explosive environment and dry place! Any other use altering from this is considered improper.

In case of damage resulting from improper use the manufacturer cannot be held liable for such risks, only the user.

Proper use also includes keeping operation, maintenance, and repair conditions specified by the manufacturer.

Maintenance and repair of the appliance must be carried out only by the manufacturer's authorized service center or his authorized service agent.

Dropping or breaking the product is not considered as intended use. Such damages caused by unintended use can not be repaired under warranty. In this case, please contact the manufacturer directly for further technical support. Damages caused by unauthorized modifications to this device results in the cessation of manufacturer's warranty.

### Safety Instructions

It is important to always keep the following safety regulations.

POLARISOL is not a medical device and it is not intended for preventing, diagnosing, monitoring, treating a disease or alleviating the symptoms of a disease. It is not suitable for the diagnosis, monitoring, treatment or alleviation of injury or handicap. It does not investigate, replace, compensate or modify the anatomy or a physiological process. POLARISOL has been exclusively developed for regeneration, well-being and cosmetic purposes.

During polarized light facial treatment, the eye area does not require any special protection. Please keep eyes always closed during treatment and if you are using contact lenses, remove them!

If you encounter any operational problems (e.g: brightness fluctuations (flashing lights) immediately pull the plug from the socket! It is important that the cable and connector will always be readily available!

Never install a damaged device (broken housing, damaged cable or optical unit)! If any liquid or other particles enter the device, immediately unplug it and contact our authorized service center!

Do not use the device in a humid or in a warmer than 30°C or colder than 10°C environment! Protect your device from scratches, dust, vibration, shock, falling!

Do not cover or stick down the ventilation slits and pay attention to that they are never blocked!

The paper packaging of the device is recyclable, it must be disposed of according to local regulations. In case the device is used properly in accordance with the specifications it requires no maintenance on your part.

The manufacturer recommends that in order to clean the internal structure of the device visit the authorized service center in every second year.

## MAINTENANCE

Appliances manufactured by Black Care Ltd. can be serviced, maintained, repaired, or any of its spare parts replaced only by the manufacturer's authorized service center or his authorized service agent.

Spare parts can only be purchased from the manufacturer.

### Other safety precautions

UNPLUG THE DEVICE BEFORE MAINTENANCE!

MAINTENANCE OF THE APPLIANCE CAN BE DONE ONLY BY THE MANUFACTURER'S AUTHORIZED SERVICE CENTER OR HIS AUTHORIZED SERVICE AGENT.

IN CASE THE POWER CORD IS DAMAGED, ELIMINATE HAZARDS BY MAKING REPLACEMENTS ONLY BY THE MANUFACTURER OR BY HIS AUTHORIZED SERVICE CENTER.

THIS APPLIANCE IS NOT INTENDED FOR USE BY PERSONS (INCLUDING CHILDREN) WITH REDUCED PHYSICAL, SENSORY OR MENTAL CAPABILITIES, OR LACK OF EXPERIENCE AND KNOWLEDGE, UNLESS THEY HAVE BEEN GIVEN SUPERVISION OR INSTRUCTION CONCERNING USE OF THE APPLIANCE BY A PERSON RESPONSIBLE FOR THEIR SAFETY.  
CHILDREN SHOULD BE SUPERVISED TO ENSURE THAT THEY DO NOT PLAY WITH THE APPLIANCE.

### Product Technical Data

	Treatment area (cm <sup>2</sup> )	Power supply	Power Consumption
POLARISOL-Pro	402	100-240 V, 50/60 Hz	24V 70W
POLARISOL-Mobilo	43	100-230 V, 50/60 Hz	12V 20W

	Wavelength (nm)	Degree of polarization
POLARISOL-Pro	400-850	>99%
POLARISOL-Mobilo	400-850	>99%

### Technical specifications POLARISOL Pro

Manufacturer:	Black Care Kft. 1139 Budapest, Rozsnyai u.11.
Model:	POLARISOL Pro
Length:	278 mm
Width:	228 mm
Height:	37 mm
Weight:	1,2 kg
Rated power:	70 W
Number of white LED:	16 pcs CREE LED
Number of infra LED:	9 pcs CREE LED
Power Supply Unit:	MEAN WELL GS90A24-P1M, 24V, 3,75A, 90W max.
Protection Class:	IP20
Power supply:	230 V 50 Hz
Rated power:	70VA 
Védettség:	IP20

POLARISOL is not a medical device and it is not intended for preventing, diagnosing, monitoring, treating a disease or alleviating the symptoms of a disease. It is not suitable for the diagnosis, monitoring, treatment or alleviation of injury or handicap. It does not investigate, replace, compensate or modify the anatomy or a physiological process. POLARISOL has been exclusively developed for regeneration, well-being and cosmetic purposes.

### Product Technical Data

Total rated power of LED bulbs	Protection class (Class II, IP20)	Weight (kg)	Operation, storage temperature	
			Operation temperature (°C)	Storage temperature (°C)
54	Class II, IP20	1,2 kg	between +10 and +30°C	between +5 and +40 °C
11	Class II, IP20	0,3 kg	between +10 and +30°C	between +5 and +40 °C

Average energy density (mW/cm <sup>2</sup> )	Light energy density per minute (J/cm <sup>2</sup> )
40	2,4
40	2,4

### Technical specifications POLARISOL Mobilo

Manufacturer:	Black Care Kft. 1139 Budapest, Rozsnyai u.11.
Model:	POLARISOL Mobilo
Length:	148 mm
Width:	88 mm
Height:	25 mm
Weight:	0,3 kg
Rated power:	20 W
Number of white LED:	4 pcs CREE LED
Number of infra LED:	2 pcs CREE LED
Power Supply Unit:	MEAN WELL GS25A12-P1J; 12V, 2,08A, 25W max.
Protection Class:	IP20
Power supply:	230 V 50 Hz
Rated power:	20VA 
Védettség:	IP20

POLARISOL is not a medical device and it is not intended for preventing, diagnosing, monitoring, treating a disease or alleviating the symptoms of a disease. It is not suitable for the diagnosis, monitoring, treatment or alleviation of injury or handicap. It does not investigate, replace, compensate or modify the anatomy or a physiological process. POLARISOL has been exclusively developed for regeneration, well-being and cosmetic purposes.

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**Customer Service:**

H - 1095 Budapest,  
Máriássy u. 7. (Ground floor)  
Tel.: +36-30-399-1099

As for the the EU Member States customers the respective law (within the meaning of the EC. 2002/96 - Directive of used electronic and electrical devices – and the validating national legislation) forbids you as the owner of this product to dispose electrical or electronic device or throw its electrical and electronic accessory elements unsorted into household waste. Please always use free delivery disposal to authorized sites for this purpose!

**© Copyright: All rights reserved.**

The content, design, photographs, illustrations used in this manual are the property of Black Care Ltd.

This publication in part or in whole – may not be reproduced, stored on a data record, transferred or distributed in any form or by means, electronic or mechanical design, copied, recorded or otherwise, without the prior written permission of Black Care Ltd.

POLARISOL is not a medical device and it is not intended for preventing, diagnosing, monitoring, treating a disease or alleviating the symptoms of a disease. It is not suitable for the diagnosis, monitoring, treatment or alleviation of injury or handicap. It does not investigate, replace, compensate or modify the anatomy or a physiological process. POLARISOL has been exclusively developed for regeneration, well-being and cosmetic purposes.

**POLARISOL**<sup>®</sup>  
POLARISED LIGHT TECHNOLOGY

POLARISOL Kft.